



Factory calibrated at 27" (the average length of the human arm), a precise 3" x 6" rectangular light pattern is produced. Forget adjusting the light due to the hot spots and color washing common with halogen lights - uniform intensity throughout the Helios' LED light pattern means illumination where you need it: throughout the oral cavity.

Reflective technology lets you work around your patient instead of your light. Since the light source is reflecting from the continuous curve of the reflector, the operating field is illuminated from all angles. No matter where your hand or instrument is, the field is being illuminated and shadows are minimized.

RANGE OF MOTION




Horizontal



Vertical



3D

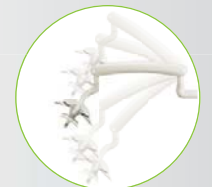

Vertical Adjustment



CONFIGURATIONS

Designed to accommodate any dental operator, the **HELIOS 3000 & 1800** are available in seven mounting variations

- 1**
CABINET
- 2**
CEILING
- 3**
WALL
- 4**
ELLIPSE



5

TRACK

6

PMU

7

POST MOUNT



(with Flexible Monitor Mount)

EFFICIENT AND RELIABLE

The Helios 3000 & 1800 consume about 70% less energy than standard halogen operating lights. Coupled with a minimum life expectancy of 30,000 hours, feel confident knowing that you are saving valuable energy and that you will never need to replace another halogen bulb again.

COOL OPERATION

Helios 3000 & 1800 are always cool to the touch, thanks to LED technology. Inherently cooler than halogen, the LEDs used in Helios mean no risk of burns for you and your staff, and no uncomfortable radiant heat.