

# Accutorr<sup>®</sup> 7

VITAL SIGNS MONITOR

## Features and Benefits

- Large 8.4" back-lit high resolution TFT display provides clear view
- Intuitive touchscreen offers easy access to information
- Manual, automatic and customizable NIBP measurement modes offer flexibility for various clinical applications
- Standard MEWS (Modified Early Warning Score) offers effective support for clinical decisions
- Manual input of vitals such as respiration rate, pain level, glucose, fluids and degree of consciousness
- Choice of hardwired or wireless EMR connectivity options
- Integration with ADT and results messaging easily accomplished via Mindray's powerful eGateway




## The Accutorr 7 Vital Signs Monitor

From spot check monitoring to continuous bedside monitoring, the Accutorr 7 not only measures NIBP, SpO<sub>2</sub> and temperature, but also has the capability for manual entries of up to six additional vitals, an early warning scoring system that can aid in clinical decisions and customizable NIBP modes for added flexibility. Every unit comes standard with these features as well as a recorder, an 8-hour lithium ion battery and all the accessories needed to begin monitoring.

In today's healthcare environment, vital sign data collection, validation and timely EMR posting are essential in providing the best possible medical care. When done manually, this process can be subject to errors and omissions which can affect efficiency and delay in clinical decisions. The Accutorr 7 offers multiple connectivity solutions to respond to your institutional requirements. Data can be posted directly to the EMR at the point-of-care, or upon completion of your rounds.

# Accutorr® 7

- Available SpO<sub>2</sub> technology includes Masimo SET® (standard) or Nellcor® Oximax®
- PI (Perfusion Index) from SpO<sub>2</sub> can guide caregivers to best measurement location
- Optional Mindray SmarTemp™ or Welch Allyn® SureTemp® Plus temperature
- Optional barcode scanner can speed patient ID input
- Stores up to 5000 measurements, accessible by patient ID for easy review
- Supports nurse call interface
- Standard RS-232 interface for third party connectivity solutions

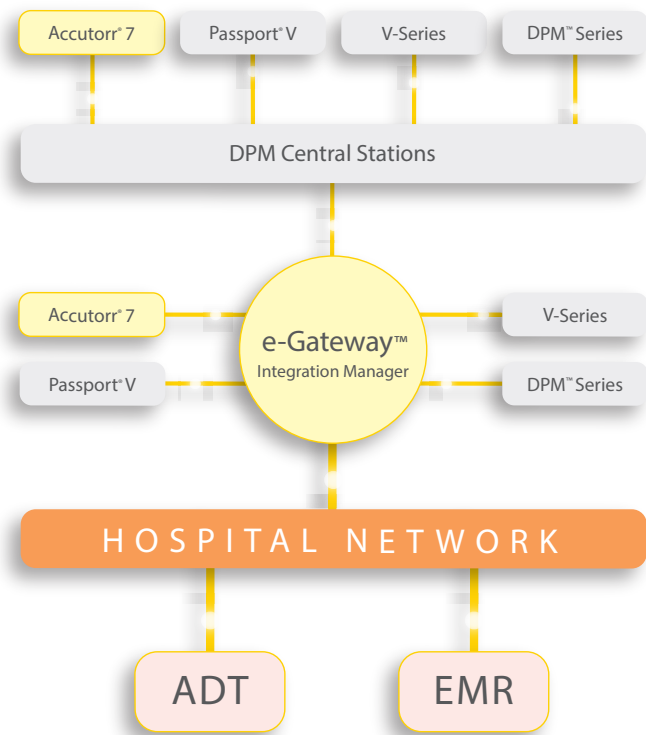

## Designed to be Connected

The Accutorr 7 can be connected to your EMR via the DPM Central Station or directly via the eGateway

The 'Manual Inputs' window contains the following data:

Pain Level	8	Temp	103.0	°F
LOC	Confused	I/O Fluid	1200 / 1000	ml
RESP	20	Glucose	2.22	mg/dl
Modifier				
Patient Position	Lying	Temp Position	Oral	
NIBP Location	L-Arm			

Buttons: Ok, Cancel

## Manual Inputs

Additional vitals and modifiers can easily be entered, saved, and sent to the EMR

The 'MEWS Scoring' window displays the following data:

PR	80	0
RESP	20	1
NIBP(mmHg)	166 / 100	0
Temp(°F)	100.4	0
AVPU	Unresponsive	3

Buttons: Calculate, Clear, MEWS Protocol >>, Record, Review

MEWS Score: 4

## Clinical Decision Support

Standard MEWS (Modified Early Warning Score) offers support for clinical decision making

Mindray DS USA, Inc.  
 800 MacArthur Blvd., Mahwah, NJ 07430  
 Tel: 800.288.2121 Tel: 201.995.8000 Fax: 800.926.4275 [www.mindray.com](http://www.mindray.com)

©2015 Mindray DS USA, Inc. Mindray® is a registered trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. All brands and product names are trademarks of their respective owners. Specifications subject to change without notice. Contact your Mindray Representative for the most current information. 05/15 P/N: 0002-08-6373 Rev C

